Silk'n®

Beautiful Technology

LED NECK MASK

ND-FM03

EN USER MANUAL

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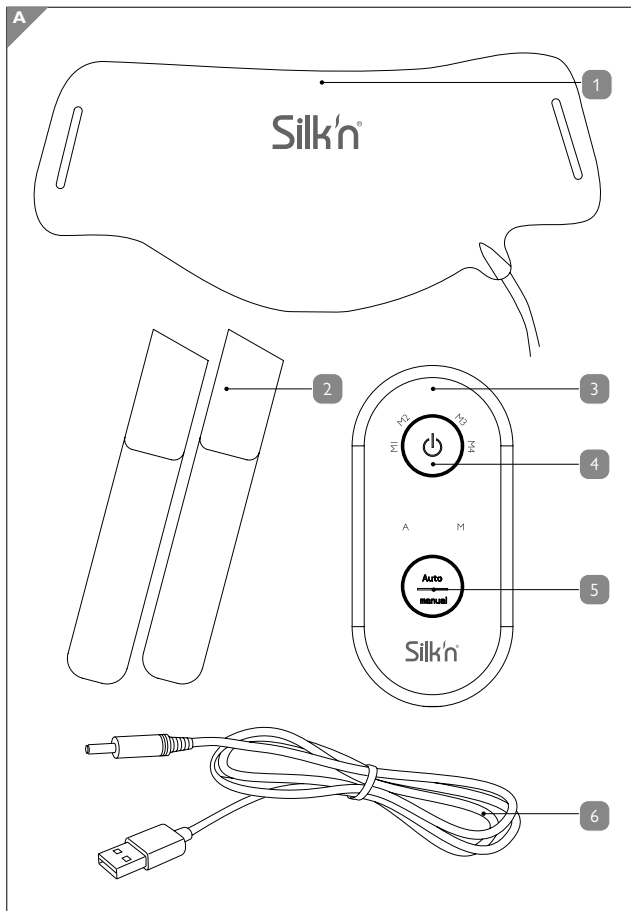
EN

DE

FR

NL

ES

Package contents / Lieferumfang / Contenu de la livraison / Leveringspakket / Volumen de suministro

Legal notice

Invention Works B.V. reserves the right to make changes to its products or specifications to improve performance, reliability, or manufacturability. Information furnished by Invention Works B.V. is believed to be accurate and reliable at the time of publication. However, Invention Works B.V. assumes no responsibility for its use. No license is granted by its implication or otherwise under any patent or patent rights of Invention Works B.V.

No part of this document may be reproduced or transmitted in any form or by any means, electronic or mechanical, for any purpose, without the express written permission of Invention Works B.V. Data is subject to change without notification.

Silk'n Beauty Ltd. has patents and pending patent applications, trademarks, copyrights, or other intellectual property rights covering subject matter in this document. The furnishing of this document does not give you any license to these patents, trademarks, copyrights, or other intellectual property rights except as expressly provided in any written agreement from Silk'n Beauty Ltd. Specifications are subject to change without notice.

Silk'n and the Silk'n logo are registered trademarks of Silk'n Beauty Ltd.

Invention Works B.V.

Galileïstraat 17, 3029 AL Rotterdam, The Netherlands

www.silkn.eu – info@silkn.eu

Device parts

1. LED Neck Mask
2. Elastic straps (2x)
3. Remote control
4. Power button
5. Auto/Manual button
6. USB cable

General information

These operating instructions accompany this LED Neck Mask ND-FM03 (referred to below only as the “product”). It contains important information about start-up and handling. Before using the product, please read through the operating instructions carefully. This particularly applies to the safety instructions. You can also find these operating instructions in PDF format on our website www.silkn.eu.

Explanation of symbols

The following symbols and signal words are used in this user manual, on the product or packaging.

WARNING!

This signal symbol/word designates a hazard with moderate risk, which may result in death or severe injury if not avoided.

NOTICE!

This signal word warns of possible damage to property.

This signal symbol provides you with useful additional information on handling and use.

Read the user manual.

Declaration of Conformity: Products labelled with this symbol meet all applicable provisions of the European Economic Area.

Devices marked with this symbol are only suitable for household use (dry indoor rooms).

This symbol identifies voltage and the current rating for direct current.

The word “GS” stands for tested safety. Products marked with this symbol meet the requirements of the German Product Safety Act (ProdSG).

Intended use

This product is intended to recover your skin's glow and elasticity. It also enhances skin metabolism, promotes anti-ageing effects, reduces inflammation and battles inflammation, acne, age spots and wrinkles. The product can be used on oily, dull and inflamed skin, but is also suitable for skin with pigmentation. This device is intended for household use only, not for commercial applications.

Safety

Only use the product as described in these operating instructions. Any other use is considered improper and may result in damage to property or persons. The manufacturer or vendor cannot be held liable for damages or injury incurred through improper or incorrect use. Before first use, please check if the product is complete. If there is any omission or if the product is not working properly, please contact customer service.

Safety instructions

- This device, including all its parts, must not be used by the following people: children, adults with limited physical, sensory or psychological capacities, those lacking experience or knowledge in how to use the device in a safe way, or those who do not understand the hazards involved.
- Children must not play with the product. Cleaning and user maintenance must not be performed by unsupervised children.
- Before treating an area for the first time, perform a skin sensitivity test by placing the

product on one of your arms. If you have no adverse reaction, start your treatment.

- Only use the included USB cable to charge the remote control. It cannot be used for other purposes. Connect the USB cable to a compatible USB power source that is easily accessible so that you can quickly disconnect it if something goes wrong.
- If the product is damaged in any way, stop using it. Contact our customer service for replacement or repair to avoid risk. Liability and warranty claims are waived in the event of repairs performed by the user, improper connection or incorrect operation.
- Do not use the product while charging.
- Keep the product and USB cable away from water or any other liquids. They are not waterproof.
- Never touch the product or USB cable with wet or damp hands.
- Stop using the product immediately if you feel any discomfort, develop any abnormal reactions (such as itching or uncomfortable inflammation) or in case your skin becomes warm to the point of discomfort.
- This product cannot be used by pregnant women, children, people with epilepsy, thyroid patients, people using thyroid medications, people who are photoallergic or use

light allergy medication.

- This product cannot be used on skin diseases, wound ulceration, trauma and post-operative wounds.
- Don't use the product for more than 15 minutes per treatment.
- **Consult your physician if you are uncertain whether or not you can use the product.**

Charging

When the battery is low, the indicator light on the remote control will flash and the product will automatically shut down after flashing 10 times.

1. Connect the remote control **3** to the mask **1**.
2. Plug the USB cable **6** into the charging port of the remote control.
3. Plug the other end of the USB cable **6** into a compatible USB power source. One to three indicator light(s) on the remote control **3** will light up. When charging is complete, all four indicator lights are lit.

Assembling

1. Connect the remote control **3** to the mask **1**.
2. Fasten the elastic straps **2** to the mask **1**.

Operation

1. Assemble the product as described in Chapter "Assembling".
2. Thoroughly clean and dry the neck and hands. Make sure they are free from any dirt.
3. Sit or lie down and place the mask **1** on your neck. Adjust the elastic straps **2** to ensure the mask fits comfortably.
4. Turn on the device and select the desired mode. The product has two modes. In the Manual Mode a preferred treatment can be selected: Red, Blue, Yellow or Purple LED light therapy. One treatment takes 10 minutes. In the Automatic Mode a preferred treatment cannot be selected. The product is activated in the following order: Red (5 min.), Blue (3 min.), Yellow (4 min.) and Purple (3 min.) LED light therapy.

Manual Mode:

- Take the remote control **3** and long press the power button **4** to turn on the device. It will enter into the Manual Mode with Red LED light therapy by default.
- Short press the power button **4** again to switch to the Blue LED light therapy.
- Short press the power button **4** a third time to switch to the Yellow LED light therapy.
- Short press the power button **4** a fourth time to switch to the Purple LED light therapy.
- Short press the power button **4** a fifth time to turn off the product. In addition, if you want to turn off the product during your treatment, you can also long press the power button **4**.
- Please note, when selecting a treatment the indicator light below the "M" flashes 5 times. This means that the product will automatically start a timer of 10 minutes. After 10 minutes the product will automatically shut down, or long press the power button **4** to shut the product off during use.

Automatic Mode:

- Take the remote control **3** and long press the power button **4** to turn on the device. It will enter into the Manual Mode with Red LED light therapy by default.
- Short press the Auto/Manual button **5** to switch to the Automatic Mode.
- Please note, when selecting the Automatic Mode the indicator light below the "A" flashes 5 times. This means that the product will automatically start a timer of 15 minutes. After 15 minutes the product will automatically shut down, or long press the power button **4** to shut the product off during use.

Don't use the product for more than 15 minutes per treatment.

5. Apply a serum and/or day/night cream after treatment.

Treatment programme

If you are using the product for the first time, please stick to the following schedule:

- First week: use once every 3 days, max. 10 minutes at a time.

- Second week: use once every 2 days, max. 10 minutes at a time.
- Third week: use once every 2 days, max. 10-15 minutes at a time.
- Fourth week: use once every day, max. 10-15 minutes at a time.

Treat your skin once or twice per week as soon as the results show.

Cleaning

Make sure to unplug the product before you clean it. Then, clean the silicone side of the mask with a water based wipe or damp clean cloth. Make sure to dry it afterwards with a dry cloth.

Storage

After cleaning, we recommend you to store the product in its original box. Store the product in a cool dry place.

Technical data

Model number:	ND-FM03
Material:	Food grade silicon & ABS
Rated voltage:	3.7V
Battery capacity:	1500mAh
Unit weight:	169 g
Number of LED lights:	57
Dimension (L × W × H):	237.5 × 108 × 8.1mm

Disposal

Disposing of the packaging

Sort the packaging before you dispose of it. Dispose of paperboard and cardboard with the recycled paper service and wrappings with the appropriate collection service.

Disposing of the product

(Applicable in the European Union and other European countries with separate collection systems of recyclable materials)

Old appliances may not be disposed of in the household waste!

Should the product no longer be capable of being used at some point in time, **dispose of it in accordance with the regulations in force in your city or county.** This ensures that old appliances are recycled in

a professional manner and also rules out negative consequences for the environment. For this reason, electrical equipment is marked with the symbol shown here.

Batteries and rechargeable batteries may not be disposed of with household waste!

As the end user you are required by law to bring all batteries and rechargeable batteries, regardless whether they contain harmful substances* or not, to a collection point run by the communal authority or to a retailer, so that they can be disposed of in an environmentally friendly manner.

Bring the entire product (including the rechargeable battery) to your collection point and make sure that it is in an uncharged state!

* labelled with: Cd = cadmium, Hg = mercury, Pb = lead

Warranty

This product is covered by a 2-year warranty according to European regulations and law. The extent of the warranty on this product is limited to technical defects caused by faulty production processes. In the event you want to claim warranty please be sure to contact our customer service for instructions. They may be able to solve your problem without needing to return the product to the store or our service centre. Our customer service is always pleased to help you!

Customer service

For more information about the products of Silk'n visit your regional Silk'n website: www.silkn.eu. If the device is damaged, defective, in need of repair or you need our assistance, please contact your nearest Silk'n service centre. This manual is also available as a PDF download from: www.silkn.eu.

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© 2022 Silk'n Beauty Ltd. All rights reserved. Silk'n and the Silk'n logo are registered trademarks of Silk'n Beauty Ltd.