PB115789A

HomeSkinovations™

www.silk'n.eu • info@silk'n.eu
InnoEssentials International
Donk 1b, 2991 LE Barendrecht
The Netherlands

Silk'n™
Beautiful Technology

Jewel

H3210

H3211

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Table of contents

1. Contraindications & warnings	24
2. Safety with Silk'n Jewel	27
3. Do not use Silk'n Jewel if:	27
4. Consult your physician before use if you have:	27
5. Get to know your Silk'n Jewel	28
6. Safety tips	29
7. Possible side effects.....	29
8. Long-term hair removal the Silk'n way	31
9. Get started.....	34
10. After treatment	36
11. Maintenance.....	37
12. Troubleshooting	37
13. Customer service.....	38
14. Warranty	38
15. Disposal	38
16. Specifications	39
17. Signs & symbols	40

1. Contraindications & warnings

Silk'n Jewel is not designed for everyone.

Hair removal by lasers or intense pulsed light sources can cause increased hair growth in some individuals. Based upon currently available data, the highest risk groups for this response are females of Mediterranean, Middle Eastern and South Asian heritage when used on the face and neck. For more information visit www.silkn.eu, www.silkn.com, or your local Silk'n domain.

- This device must not be used by the following people: children, adults with limited physical, sensory or psychological capacities, or those lacking experience or knowledge in how to use the device. Children should be supervised to ensure that they do not play with the device.
- Do not expose treated areas to the sun, and do not use energy levels 2-5 on tanned skin or after recent sun exposure (including artificial tanning machines). Refer to section 6.1 Exposure to the sun.
- Possible hazardous optical radiation is emitted from this device. Do not look at the operating lamp. Eye injury may result.

- ⦿ Use below the cheekbones and not around the eyes. Do not use on areas with injected fills or toxins.
- ⦿ Do not treat on tattoos or permanent makeup, dark brown or black spots, (such as large freckles, moles, or warts), nipples, genitals, around the anus, or lips.
- ⦿ Do not treat on active eczema, psoriasis, lesions, open wounds or infections (cold sores), abnormal skin conditions caused by systemic or metabolic diseases (diabetes for example).
- ⦿ Do not treat on abnormal skin conditions caused by diabetes, systemic or other metabolic diseases.
- ⦿ Do not treat an area that has a history of herpes outbreaks.
- ⦿ Do not use if you are pregnant or nursing.
- ⦿ Do not treat on areas where you may want hair to grow. The results are permanent.
- ⦿ Do not treat the same area of skin more than once per hair removal session.
- ⦿ Stop use immediately if your skin is very red, or if your skin blisters or burns.
- ⦿ Use a lower energy level if you experience any discomfort.

- Never use flammable liquids such as alcohol (including perfumes, sanitizers, or other applications containing alcohol) or acetone to clean the skin before use.
- Keep your device away from water or liquids! Do not use near bathtubs, showers, basins or other vessels containing water. Do not place or store the device where it can fall, be pushed or placed anywhere near water. There is a risk of electrocution. Unplug immediately if the device comes into contact with any fluids and contact Customer Service.
- Remove all jewelry before use.
- Never leave the device unattended when plugged into an outlet. Keep the AC/DC adapter accessible at all times, so it can be removed from the socket easily in case of emergency.
- Do not attempt to open or repair your Silk'n device. You may expose yourself to dangerous electrical components or pulsed light energy and cause serious harm to yourself. It will also void your warranty.

2. Safety with Silk'n Jewel

With Silk'n Jewel safety comes first

Home Pulsed Light™ (HPL) technology in Silk'n Jewel provides superior safety with lower energy levels. HPL achieves long-term hair removal results at a fraction of the energy level used in other light-based hair removal equipment. The low energy reduces the risk of harm or complications, and contributes to your overall safety.

3. Do not use Silk'n Jewel if:

- The device has a damaged cord or plug. Keep the power cord away from heated surfaces.
- You see or smell smoke when it is in use.
- The **skin colour sensor, treatment surface, or device casing** is cracked or coming apart.
- The device is not working or appears to be damaged in any other way.
- The plug, attachments or accessories for this device are not recommended by Home Skinovations Ltd.

In such cases, contact Silk'n Customer Service.

4. Consult your physician before use if you have:

- A history of keloidal scar formation.
- Epilepsy.
- A disease related to photosensitivity, such as porphyria, polymorphic light eruption, solar urticaria, lupus, etc.
- A history of skin cancer or areas of potential skin malignancies.
- Received radiation therapy or chemotherapy treatments within the past 3 months.
- Used physician-prescribed drugs on your skin within the past 6 months.

- Any other condition which, in your physician's opinion, would make it unsafe for you to be treated.

5. Get to know your Silk'n Jewel

Silk'n Jewel is a light-based device for permanent hair reduction designed for home use. It is a powerful electrical device and should be used with special attention to safety. Be sure to read **all** warnings, contraindications and safety precautions before use, and follow them strictly. We recommend you re-familiarize yourself with the instructions before each treatment.

5.1 Package contents

Your Silk'n Jewel package includes:

- Hair removal device
- Adapter with cable
- User manual
- Warranty card

5.2 Device description

Refer to the diagram at the front of this manual to see the location of each feature of your Silk'n Jewel.

1. Adapter socket
2. Energy level indicator lights (1-5)
3. Control button
4. Ready light
5. Skin colour sensor light
6. Treatment surface
7. Skin colour sensor

5.3 Intended use

Silk'n Jewel is intended to remove unwanted hair and effect long-term, or permanent hair reduction. Permanent hair reduction is defined as a long-term stable reduction in the number of hairs re-growing after a treatment regime. You may use Silk'n Jewel to remove unwanted body hair. Ideal body areas include the underarms, bikini line, arms, legs, face, back, shoulders and chest.

Do not use for any purpose other than as described in this manual.

6. Safety tips

This section has information to help you use your device safely.

6.1 Exposure to the sun

- The lowest energy level (1) is designed for safe use even after sun exposure.
- If you use energy levels 2-5, avoid further exposure to the sun for **4 weeks before and 2 weeks after** your treatment. This is because higher quantities of melanin expose the skin to greater risk when using any method of light-based hair removal. This applies to all skin types, even those that don't seem to tan quickly.
- We don't advise it, but if you have to go out in the sun after a treatment, be sure to apply sunscreen SPF30 or higher to the treated area. Do this for 2 weeks after a treatment.

6.2 Exposure of pulse light to the eyes

- Although you can use the device to treat facial hair, take extra care to avoid contact with the eyes, and use only below the cheekbones.
-  The device will only emit repeated pulses when it is in firm contact with your skin, and you keep moving it across your skin.

6.3 Wounds or infections

- Wait for the affected area to heal before further use.

7. Possible side effects

If you follow the instructions fully and pay attention to all warnings and contraindications, complications and side effects are uncommon. However every cosmetic procedure, including those designed for home use, involves some degree of risk. Therefore, it is important

that you understand and accept the risks and complications that may occur with pulsed light hair removal systems designed for home use.

Minor skin discomfort

Home pulsed light hair removal is generally well-tolerated. However, some users feel mild discomfort during use. This is best described as a mild stinging sensation on the treated skin areas. The stinging sensation may last during or after treatment.

This means that you should use Silk'n Jewel at a lower energy level, and if that doesn't help, you should stop using the device.

Skin redness

Your skin may become red immediately after using the device or within 24 hours. Redness generally clears up within 24 hours.

See your doctor if redness does not go away within 2 to 3 days.

Make sure you use the energy level that is comfortable for you and suitable for your skin colour.

Increased sensitivity of the skin

The skin of the treated area is more sensitive so you may encounter dryness or flaking of the skin.

Pigment changes

Silk'n Jewel targets the hair shaft, in particular the pigmented cells in the hair follicle and the hair follicle itself. Nevertheless there is a risk of temporary hyperpigmentation (increased pigment or brown discolouration) or hypopigmentation (whitening) to the surrounding skin. Risk of changes in skin pigmentation is higher for people with darker skin tones or heavily tanned skin. Usually discolouration or changes to skin pigment are temporary, and permanent hyperpigmentation or hypo-pigmentation are rare.

Excessive redness and swelling

In rare cases, treated skin may become very red and swollen. This is more common in sensitive areas of the body. The redness and swelling should subside within 2 to 7 days and should be treated with frequent applications of ice. Gentle cleansing is OK, but avoid exposure to the sun.

Blisters or burns

Cool the affected area and apply anti-burn cream. Wait until the area has healed completely before another treatment, and use a lower energy level. **If the skin does not return to normal within 3 days, consult your physician.**

8. Long-term hair removal the Silk'n way

The efficiency of hair removal varies from person to person according to body area, hair colour, and how Silk'n Jewel is used. A typical full hair growth cycle may take 18-24 months. During this time, multiple treatment sessions are required to achieve permanent hair removal.

A typical Silk'n Jewel hair removal plan during a full hair growth cycle includes:

- Treatments 1-4: Plan two weeks apart.
- Treatments 5-7: Plan four weeks apart.
- Treatments 8+: Treat as needed, until desired results are achieved.

8.1 Expected treatment outcomes

Individual response depends on hair type as well as biological factors that may affect hair growth patterns. Some users may respond quicker or slower than the average number of treatments.

Skin Type ¹	Body Hair ²	Light Fluence [J/cm ²]	Part of the Body	Average Treatment Time	Effectiveness of Hair Removal
1-5	Dark blonde to brown	2.5-4	Lower or upper legs		60%
			Arms	Session 10-12, approx. 4 weeks apart.	50%
			Bikini line		60%
			Underarms		55%

Skin Type ¹	Body Hair ²	Light Fluence [J/cm ²]	Part of the Body	Average Treatment Time	Effectiveness of Hair Removal
1-5	Brown to black	2.5-4	Lower or upper legs	Sessions 8-10 approx. 4 weeks apart.	70%
			Arms		60%
			Bikini line		70%
			Underarms		65%

¹ Skin types are based on the Fitzpatrick scale of skin colours. Refer to the skin colour chart in the diagram at the front of the manual.


² Not effective on white body hair.

8.2 Energy levels

The energy level determines the intensity of the light pulse delivered to your skin. As the energy level increases, so do hair removal results, as well as the risk of possible side effects and complications.

- The device has 5 levels of energy from 1 (the lowest) to 5 (the highest), represented by 5 **energy level indicator lights**.
 - When the device is turned on, the energy level is set to level 1 automatically. Only one **energy level indicator light** will be on.
 - To change the energy level, press the **control button** repeatedly to increase the energy level from 1 to 5. Once at level 5, press the **control button** again to reduce the energy level from 5 back to 1.
-  To determine the right energy level for your treatment area, refer to section 9.1 Before your first treatment.

8.3 Pulsing and gliding

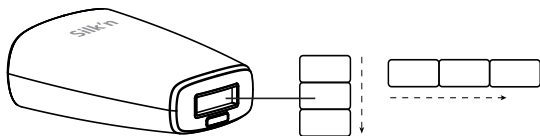Pulsing is the best technique for treating specific spots or small areas.

- Press the **treatment surface** against your skin. Work in rows (either across or down as shown in the image). Start at one end of each row and progress to the other end.
- The **treatment surface** creates temporary pressure marks on the treated area so you can see areas you have treated.
- After each pulse, move the **treatment surface** to the next spot. Keep doing this pulsing action for the area you are treating.
- Avoid gaps or overlaps as you use the device.

Gliding is the best technique for treating larger parts of the body.

- Press the **treatment surface** firmly against your skin and slowly glide over your skin (either across or down as shown in the image). Make sure that the **treatment surface** is always touching your skin.
- At lower energy levels, the device takes only 1-1.5 seconds to recharge between pulses. This generates a nonstop sequence of flashes.
- Glide at a speed so that you avoid gaps or overlaps.

⚠ Do not treat the same area of skin more than once per hair removal session.

⚠ Stop use immediately if you have any discomfort.

8.4 The skin colour sensor

Light-based hair removal on darker skin tones may result in adverse effects such as burns, blisters, and skin colour changes (hyper- or hypo-pigmentation). Your Silk'n Jewel has a **skin colour sensor** that detects the colour of your skin automatically. The skin colour sensor checks your skin colour at the start of, and occasionally during, each treatment.

-  If your skin tone is too dark for safe use, the device will not operate, and the skin colour sensor light will blink.

8.5 What to expect when treating with Silk'n

For many people, using Silk'n Jewel may be their first experience with a light-based device. During a treatment session it is normal to experience or feel:

- **Noise:** The cooling fan makes noise similar to a hair dryer. This is normal.
- **A popping sound:** With each flash of HPL™ light.
- **A sensation of warmth and tingling:** During each light pulse.
- **Some redness:** It is not uncommon to see some very mild, pink-like skin colouration during or just after a treatment. This is usually most noticeable around the hairs themselves.

 Stop use immediately if you see your skin is very red.

 Contact your physician immediately, if you experience any complication.

9. Get started

To use your Silk'n Jewel, follow these steps:

1. Remove the device and other components from box.
2. Plug the **adapter plug** into the device socket and plug the **adapter** into an electrical outlet. All **energy level indicator lights** will turn on, and then go off.

9.1 Before your first treatment

Remove all jewellery before use.

You must decide which energy level is right for your skin, and the area you want to treat. The device has 5 levels of energy from 1 (the lowest) to 5 (the highest). For example, you might treat one area at energy level 1, and another at level 4. For information on energy levels, refer to section 8.2 Energy levels.

To choose the right energy level, test each treatment area separately (refer to the skin colour chart at the front of this manual).

1. Press the **control button** to turn on the device. The level 1 **energy level indicator light** will go on.
2. Place the **treatment surface** on your skin to deliver one pulse. If you do not experience abnormal discomfort press the **control button** again. This will increase the energy level. Deliver one pulse on a different spot at energy level 2.
3. Continue testing at increased energy levels. Use the energy levels indicated on the skin colour chart. When you find the highest energy level that is comfortable for that treatment area, you can start your treatment.

Use a lower energy level if you experience any discomfort.

Repeat this test for each part of the body you wish to treat.

You are now ready to use your Silk'n Jewel

Before treating an area, the skin should be shaved, clean, dry and free of any powders, antiperspirants or deodorants. **Do not wax, pluck or tweeze the hair out.**

1. Be sure to cover up moles, warts or other marks on the skin that should not be treated.
2. Press the **control button**. The device will turn on, the level 1 **energy indicator light** will be on, and a fan sound will start.
3. Approximately 1 second after pressing the **control button**, the **ready light** will turn on. The device is then ready to use at the lowest energy level.

4. **Set the energy level.** Press the **control button** repeatedly to increase the energy level until you have the level you want. The **energy level indicator lights** indicate the energy level setting (refer to section 8.2 Energy levels).
5. Decide if you wish to treat the area by **Gliding** or **Pulsing** (refer to section 8.3 Pulsing and gliding).
6. Apply the **treatment surface** to the skin. Make sure that the skin is spread evenly and smoothly. When the **treatment surface** is in full contact with the skin, the **ready light** will start blinking slowly.
7. The device will first determine the colour of your skin automatically and it will flash a pulse of light onto your skin. You will feel a mild sensation of warmth and tingling, and you may hear a popping sound.
8. If you are pulsing, remove the device from your skin. The device will immediately recharge for the next pulse. After a few seconds (depending on the energy level) the **ready light** will turn on again, and start blinking slowly when pressed against the skin.
9. If you are pulsing, repeat the process. If you are gliding, simply continue. Continue to pulse or glide until you have treated the area fully.

Stop use immediately, if your skin blisters or burns.

10. After treatment

When you have finished your treatment session:

1. Turn the device off by giving the control button an extended press. Remember the last energy level setting you used, because it will not be restored when you turn the device on again.
2. Remove the adapter from the electrical outlet.

11. Maintenance

It is important to clean your Silk'n Jewel device after each treatment.

- Unplug the device.
 - To clean the device including the treatment surface, wipe with a clean, dry cloth. Never use acetone or other abrasive fluids, because they will damage your device.
 - After cleaning, we recommend you store the device in its original box.
-  Never immerse your Silk'n Jewel or any of its parts in water.

12. Troubleshooting

Do not attempt to open or repair your Silk'n device. You may expose yourself to dangerous electrical components or pulsed light energy and cause serious harm to yourself. It will also void your warranty.

Problem	Checks
My device does not start.	Make sure the adapter is properly connected to the device and that it is plugged into an electrical outlet on the wall.
There is no light pulse when the device touches my skin.	Check that the device is connected and the ready light is on. Make sure that the treatment surface is evenly and firmly pressed to the skin and that you are moving the device across your skin. It is also possible that your skin colour is too dark.
The skin colour sensor light and the ready light blink.	The skin colour is too dark.
One or more lights blink.	The device has malfunctioned.

Problem	Checks
The device has stopped working and is in standby mode.	Press the control button again to start the device.

For problems or queries, please contact Silk'n Customer Service.

13. Customer service

See our website for more information: www.silkn.eu. This manual is also available as a PDF download from: www.silkn.eu. Contact Customer Service to report unexpected operations or events, or any other problem with your device. Our contact details are:

Country	Service number	Email
Belgium	0900-25006	serviceconsommateurbe@silkn.eu
France	0891-655557	serviceconsommateurfr@silkn.eu
Germany	089 51 23 44 23	kundenservice@inno-essentials.de
Italy		servizioclientiit@silkn.eu
Netherlands	0900-2502217	klantenservicenl@silkn.eu
Spain	900 823 302	servicioalcliente@silkn.eu
UK	0906-2130009	customer-careuk@silkn.eu
Other countries	+31 (0)180-330 550	info@silkn.com

14. Warranty

Refer to your Silk'n warranty card.

15. Disposal

If your device reaches its End of Life stage, be sure to follow your country's laws for the disposal of electrical items.

Dispose of the packaging separated into single type materials. Dispose of cardboard and carton as waste paper and film via the recyclable material collection service.

Do not dispose with household waste! The WEEE directive requires that you dispose of this electrical device (including all accessories and parts) at a designated recycling collection point for electrical devices. This ensures that the device is recycled professionally, and prevents the release of harmful substances* into the environment. Make sure that the device is in an uncharged state before disposal.

* Products labelled with Cd = cadmium, Hg = mercury, Pb = lead.

16. Specifications

Model no.	H3210 / H3211
Spot size	0.9cm x 3cm [2.7cm ²]
Technology	Home Pulsed Light
Maximum energy level	Max 4J/cm ²
Wavelength	475-1200nm
Operation & safety	Skin colour sensor ensures use on appropriate skin types
System weight	200 gm
Transport & storage between uses and storage condition	Temperature: -40°C to +70°C Relative humidity : 10% to 90% Atmospheric pressure: 500 to 1060 hPa
Operating conditions	Temperature: 10°C to 35°C Relative humidity: 30% to 75% Atmospheric pressure: 700 to 1060 hPa
Time from minimum storage temperature between uses to minimum operating temperature with ambient temperature of 20°C degrees	30 minutes

Time from maximum storage temperature between uses to maximum operating temperature with ambient temperature of 20°C degrees	30 minutes
Mode of operation	Continuous
Adapter models and rating	KSA24A1200150Hy PSM18x-120L6HS For y, refer to the country / region
Input Output	100-240 VAC, 50-60 Hz, 0.5A 12 VDC, 1.5A
Expected service life of the equipment	5 years

17. Signs & symbols

	Products marked with this symbol fulfil the requirements of the EU directives.
	Class II. The adapter is protected against electric shock by means of double insulation.
IP30	The device is protected against the ingress of objects larger than 2.5mm diameter, such as fingers, to protect against electric shock.
	Do not use in wet environments.
	Do not dispose with household waste. Dispose in accordance with your country's laws and help protect the environment. In accordance with the WEEE Directive.
	Conforms to technical regulations required by the Eurasian Customs Union.

Legal Notice

Home Skinovations Ltd reserves the right to make changes to its products or specifications to improve performance, reliability, or manufacturability. Information furnished by Home Skinovations Ltd is believed to be accurate and reliable at the time of publication. However, Home Skinovations Ltd assumes no responsibility for its use. No license is granted by its implication or otherwise under any patent or patent rights of Home Skinovations Ltd.

No part of this document may be reproduced or transmitted in any form or by any means, electronic or mechanical, for any purpose, without the express written permission of Home Skinovations Ltd. Data is subject to change without notification.

Home Skinovations Ltd has patents and pending patent applications, trademarks, copyrights, or other intellectual property rights covering subject matter in this document. The furnishing of this document does not give you any license to these patents, trademarks, copyrights, or other intellectual property rights except as expressly provided in any written agreement from Home Skinovations Ltd. Specifications are subject to change without notice.

Silk'n and the Silk'n logo are registered trademarks of Home Skinovations Ltd, Tavor Building, Shaar Yokneam P.O.Box 533, Yokneam 2069206, ISRAEL
www.silkn.eu or www.silkn.com, info@silkn.com.

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