
beurer

HK 68 Cosy

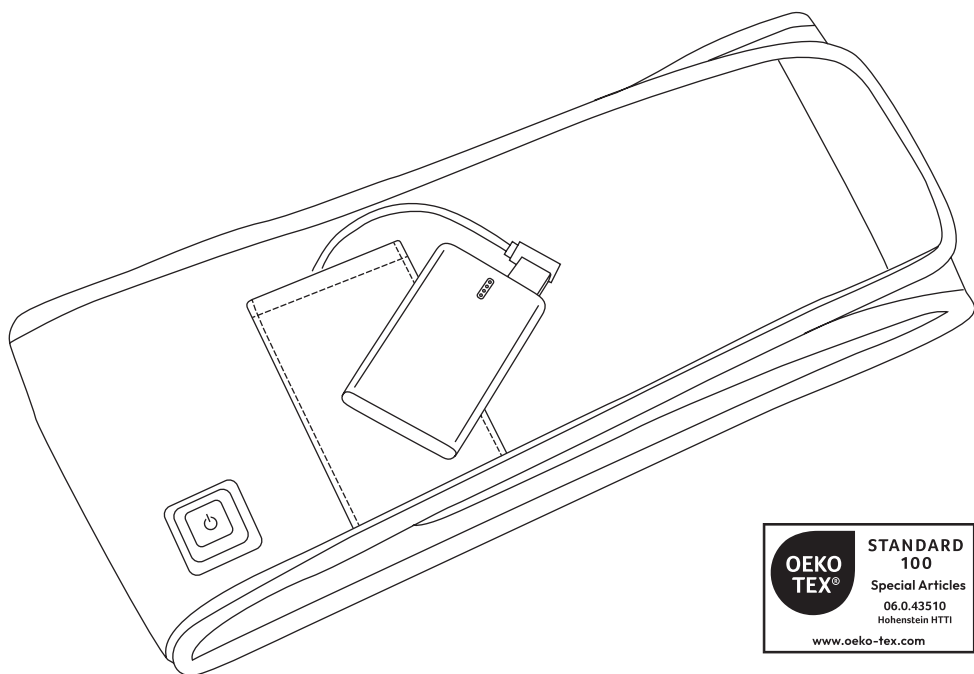DE Akku-Heizgürtel Gebrauchsanweisung	2	PL Pas rozgrzewający z akumulatorem Instrukcja obsługi	29
EN Battery powered heat belt Instruction for Use.....	9	CS Vyhřívací bederní pás s akumulátorem Návod k použití	36
FR Ceinture chauffante à batterie Mode d'emploi	15	SK Akumulátorový vyhrievací pás Návod na použitie	42
TR Akulu ısı kemer Kullanım kılavuzu	22	HU Elektromos derékmelegítő öv, akkus Használati utasítás	48

DEUTSCH

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Read these instructions for use carefully. Observe the warnings and safety notes. Keep these instructions for use for future reference. Make the instructions for use accessible to other users. If the device is passed on, provide the instructions for use to the next user as well.

Contents

1. Included in delivery	10	4.5 Switching on	13
1.1 Device description	10	4.6 Setting the temperature	13
2. Important instructions		4.7 Automatic switch-off time	13
Retain for future use	10	4.8 Switching off	13
3. Intended use	12	5. Cleaning and maintenance	13
4. Initial use / Operation	12	6. Storage	14
4.1 Safety	12	7. What if there are problems?	14
4.2 Charging the battery pack	12	8. Disposal	14
4.3 Checking the state of charge of the battery pack	12	9. Technical specifications	14
4.4 Operating the heat pad	13	10. Guarantee/service	14

Signs and symbols

The following symbols are used on the device, in these instructions for use, on the packaging and on the type plate for the device:

	WARNING: Warning of risks of injury or health hazards		
	CAUTION: Safety information about possible damage to appliance/accessories		
	NOTE: Important information		
	Read the instructions		Do not insert pins
	Do not use folded or rucked		Not to be used by very young children (0 - 3 years)
	Manufacturer		
	Protection class III device		
	Please dispose of the device in accordance with EC Directive – WEEE (Waste Electrical and Electronic Equipment).		
	Batteries must not be disposed of with household waste		
	Marking to identify the packaging material. A = Material code, B = Material number: 1-7 = Plastics, 20-22 = Paper and cardboard		
	This product satisfies the requirements of the applicable European and national directives.		
	This device meets the high requirements of OEKO-TEX® STANDARD 100 for human ecology as verified by the Hohenstein Research Institute.		


	Wash at a maximum temperature of 30 °C, Very gentle wash		Do not iron
			Do not bleach
	Do not dry-clean		Do not dry in a tumble dryer

1. INCLUDED IN DELIVERY

- 1 x Cordless heated belt
- 1 x Battery pack
- 1 x Charging cable
- 1 x Instructions for use

1.1 Device description

1. Battery pack
2. ON/OFF button of the battery pack
3. Connection socket
4. Connection socket for the battery pack
5. LED display for the charge level (white)
6. ON/OFF button of the heat belt with temperature setting display
7. Connector plug for the connection cable
8. Charging cable

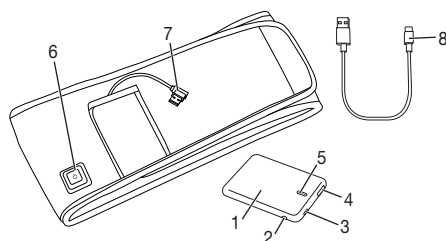2. IMPORTANT INSTRUCTIONS RETAIN FOR FUTURE USE

WARNING

- Non-observance of the following notes may result in personal injury or material damage (electric shock, skin burns, fire). The following safety and hazard information is not only intended to protect your health and the health of others, it should also protect the product. For this reason, pay attention to these safety notes and include these instructions when handing over the product to others.
- This heat pad must not be used by people who are not sensitive to heat or by vulnerable people who may not be able to react to overheating (e.g. diabetics, persons with skin changes related to illness or persons with scarred skin in the area where the pad is to be applied, after taking any pain-relieving medication or alcohol).
- This heat pad must not be used by very young children (0-3 years old) as they are unable to respond to overheating.
- The heat pad can be used by children older than 3 and younger than 8 years of age provided they are supervised. For this, the control must always be set to the minimum temperature.

- This heat pad may be used by children over the age of 8 and by people with reduced physical, sensory or mental skills or a lack of experience or knowledge, provided that they are supervised or have been instructed on how to use the heat pad safely, and are fully aware of the consequent risks of use.
- Children must not play with the heat pad.
- Cleaning and user maintenance must not be performed by children unless supervised.
- This heat pad is not designed for use in hospitals.
- The device is only intended for domestic/private use, not for commercial use.
- Do not insert pins.
- Do not use folded or rucked.
- Do not use if wet.
- This heat pad may only be used in conjunction with the battery pack specified on the label.
- If the heat pad is being used over several hours, we recommend that you set the lowest temperature setting to avoid overheating the heated body part, which may lead to burns to the skin.
- The electrical and magnetic fields emitted by this heat pad may interfere with the function of a pacemaker. However, they are still well below the limits: electrical field strength: max. 5000 V/m, magnetic field strength: max. 80 A/m, magnetic flux density: max. 0.1 millitesla. Please consult your doctor and the manufacturer of your pacemaker before using this heat pad.
- Do not pull, twist or make sharp bends in the cables.
- Please check this heat pad frequently for signs of wear and tear or damage. If any such signs are evident, if the heat pad has been used incorrectly or if it no longer heats up, it must be checked by the manufacturer before being switched on again.
- Under no circumstances should you open or repair the heat pad (including the accessories) yourself because faultless functionality can no longer be guaranteed thereafter. Failure to observe this will invalidate the guarantee.
- If the cables of the heat pad or the charger are damaged, they must be replaced by the manufacturer, its customer service department or a similarly qualified person to avoid any risk.

- Do not expose control and cables to direct sunlight. Otherwise the heat pad may get damaged.
- When this heat pad is switched on:
 - Do not place any sharp objects on it
 - Do not place any heat sources, such as hot water bottles, heat pads or similar, on it.
- While charging, do not cover the battery pack with covers, cushions or similar items.
- The cells in the rechargeable battery pack cannot be replaced.
- Always observe the notes on usage, cleaning and maintenance and storage.
- Allow the heat pad to cool down before storing it away. Otherwise the heat pad may be damaged.
- To avoid making sharp bends in the heat pad, do not place any objects on top of it whilst storing.
- If you should have any questions about using our devices, please contact Customer Services.

3. INTENDED USE

CAUTION

This heat pad is only designed to warm up the human body.

4. INITIAL USE / OPERATION

4.1 Safety

CAUTION

The heat pad is fitted with a SAFETY SWITCH-OFF. A sensor prevents the heat pad from overheating by means of automatic switch-off in the event of an error. Once the SAFETY SWITCH-OFF has taken place, the ON/OFF button on the heat pad and the LED display on the battery pack are no longer lit up when switched on. After cooling down to room temperature, the heat pad can be switched on again.

4.2 Charging the battery pack

CAUTION

The battery pack (1) must be charged prior to initial use or if the device has not been used for a long period.


- To do this, plug the supplied charging cable into a suitable charger (5V, 2A). Connect the charging cable to the connection socket of the battery pack (3). Now plug the charger into a socket with the mains voltage corresponding to the specifications on the charger. During the charging process, the white LEDs on the battery pack flash according to the charging progress.
- The battery pack is fully charged when all four white LEDs are permanently lit (it takes approx. 3.5 hours to fully charge a flat battery pack).
- After the battery pack has been disconnected from the charging cable, the white LEDs light up for a few seconds depending on the charge level.
- The battery pack should be charged at least 1x per year to maintain its capacity.
- While charging, do not cover the battery pack with covers, cushions or similar items.
- Only charge the battery pack in a dry environment at normal room temperature.

4.3 Checking the state of charge of the battery pack

- Pressing the ON/OFF button on the battery pack once indicates the current state of charge of the battery pack (4 white LEDs).

4.4 Operating the heat pad

- Place the charged battery pack into the built-in pocket on the heat pad. Then insert the connector plug for the connection cable (7) of the heat pad into the connection socket on the battery pack (4). The heat pad is now ready for use.

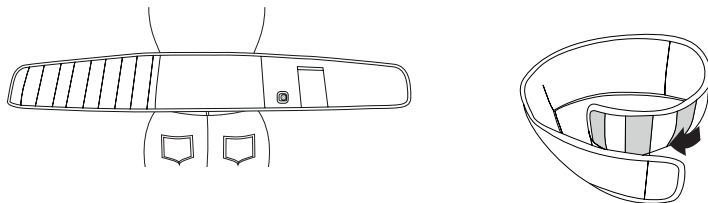- Place the heat pad on your stomach or back so that the battery pack faces outwards. Subsequently guide the right strap of the heat pad around the body and secure the heat pad by closing the hook-and-loop fastener.

IMPORTANT

Ensure that the heat pad will not bunch up or become folded during use.

4.5 Switching on

Press the ON/OFF button (6) until it lights up. It is illuminated when switched on.

4.6 Setting the temperature

To change the temperature, press the ON/OFF button (6) repeatedly.

Setting 1: minimum heat (white)

Setting 2: medium heat (blue)

Setting 3: maximum heat (red)

NOTE

The fastest way to warm up the heat pad is to initially set the highest temperature setting.

WARNING

If the heat pad is being used over several hours, we recommend that you set the lowest temperature setting to avoid overheating the heated body part, which may lead to burns to the skin.

NOTE

When the battery pack (1) has been almost completely run down, the heat pad is switched off automatically to protect the battery pack. The battery pack must then be charged.

4.7 Automatic switch-off time

This heat pad is equipped with an automatic switch-off function. This turns off the heat supply approx. 90 minutes after the initial use. The ON/OFF button on the heat pad will then no longer light up. You can switch the heat pad on again immediately by pressing and holding the ON/OFF button on the heat pad. If the heat pad can no longer be switched on, then the battery pack (1) has been run down.

4.8 Switching off

Switch off the heat pad by pressing and holding the ON/OFF button on the heat pad (6).

NOTE

Please note that the battery pack must be protected against extreme cold (such as being placed under clothing) to maintain its performance.

5. CLEANING AND MAINTENANCE

CAUTION

Prior to cleaning, you must always disconnect the battery pack (1) from the heat pad and remove it from the heat pad pouch, otherwise there is a risk of damage.

CAUTION

The battery pack (1) must never come into contact with water or other liquids, as this may cause damage.

- To clean the battery pack, use a dry, lint-free cloth. Do not use any chemical or abrasive cleaning agents.
- Small marks on the heat pad can be removed with a damp cloth and if necessary, with a little liquid detergent for delicate laundry.

CAUTION

Please note that the heat pad may not be dry cleaned, wrung out, tumble dried, put through a mangle or ironed. Otherwise the heat pad may be damaged.

- This heat pad is machine-washable. For this purpose, close the hook-and-loop fastener of the heat pad and place it in a pillow case.
- Set the washing machine to a particularly gentle wash cycle at 30 °C (wool cycle). Use a delicate laundry detergent and measure it out according to the manufacturer's instructions.
- Let the heat pad dry flat on a clothes horse.

CAUTION

- Do not use pegs or similar items to attach the heat pad to the clothes horse, otherwise the heat pad may be damaged.
- Only reconnect the heat pad to the battery pack (1) once the heat pad and the connector plug of the connection cable (7) are completely dry, otherwise this could damage them.

6. STORAGE

If you do not plan to use the heat pad for a long period of time, we recommend that you store it in the original packaging.

CAUTION

- Please allow the heat pad to cool down before storing. Otherwise the heat pad may be damaged.
- To avoid making sharp bends in the heat pad, do not place any objects on top of it whilst storing.
- The battery pack should be charged 1x per year to maintain its capacity.

7. WHAT IF THERE ARE PROBLEMS?

Problem	Cause	Solution
The ON/OFF button on the heat pad does not light up after being pressed and held when: – the battery pack is connected to the heat pad.	The battery pack has been fully run down.	Charge the battery pack.
	There was a SAFETY SWITCH-OFF due to an error.	Allow the heat pad to cool down to room temperature and switch on again.
The ON/OFF button on the heat pad does not light up after being pressed and held when: – the fully charged battery pack is connected to the heat pad.	The heat pad or the battery pack is faulty.	Contact Customer Service
The white LED display for the state of charge lights up or flashes when the battery pack is not connected to the heat pad.	Display for the state of charge of the battery pack: LED lights up for approx. 35 seconds: state of charge > 25% LED flashes for approx. 35 seconds: state of charge < 25%	Contact Customer Service

8. DISPOSAL

For environmental reasons, do not dispose of the device in the household waste at the end of its useful life. Dispose of the unit at a suitable local collection or recycling point. Dispose of the device in accordance with EC Directive – WEEE (Waste Electrical and Electronic Equipment). If you have any questions, please contact the local authorities responsible for waste disposal.

9. TECHNICAL SPECIFICATIONS

See the type plate label on the heat pad, on the battery pack and on the charger.
 Operating temperature: 0 to +45°C
 Storage temperature: -20°C to +45°C
 Recommended storage temperature: -20°C to +25°C
 Subject to technical changes

10. GUARANTEE/SERVICE

Further information on the warranty and warranty conditions can be found in the warranty leaflet supplied.

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