
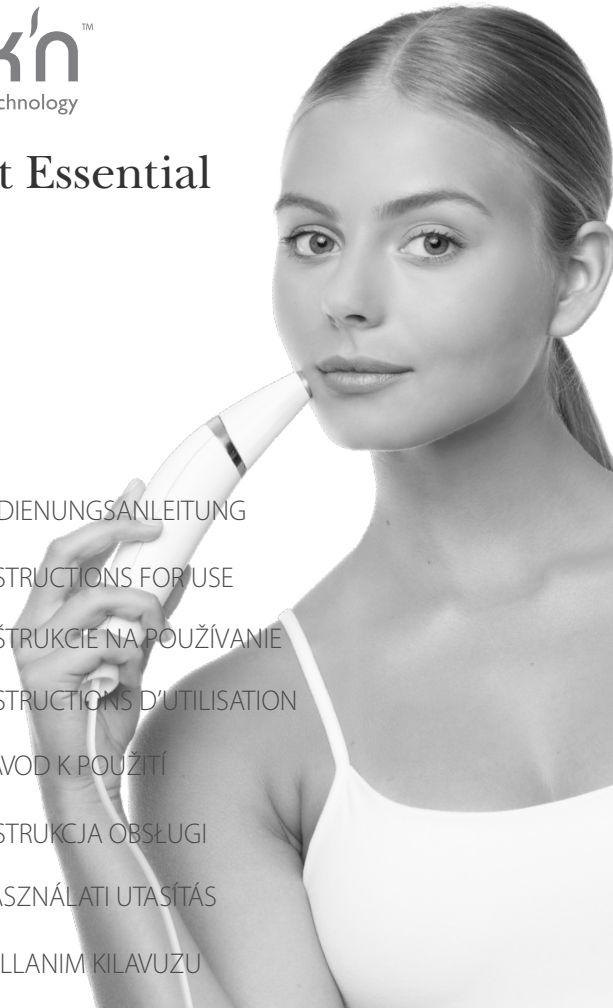
ReVit Essential

H4201

- 
- A black and white photograph of a woman with her hair pulled back, holding a white Silk'n ReVit Essential device to her cheek. The device is a handheld, ergonomic tool with a white body and a silver-colored tip. A thin white cord is attached to the bottom. The woman is looking directly at the camera with a neutral expression.
- DE BEDIENUNGSANLEITUNG
 - EN INSTRUCTIONS FOR USE
 - SK INŠTRUKCIE NA POUŽÍVANIE
 - FR INSTRUCTIONS D'UTILISATION
 - CS NÁVOD K POUŽITÍ
 - PL INSTRUKCJA OBSŁUGI
 - HU HASZNÁLATI UTASÍTÁS
 - TR KULLANIM KILAVUZU

Silk'n™

ReVit Essential

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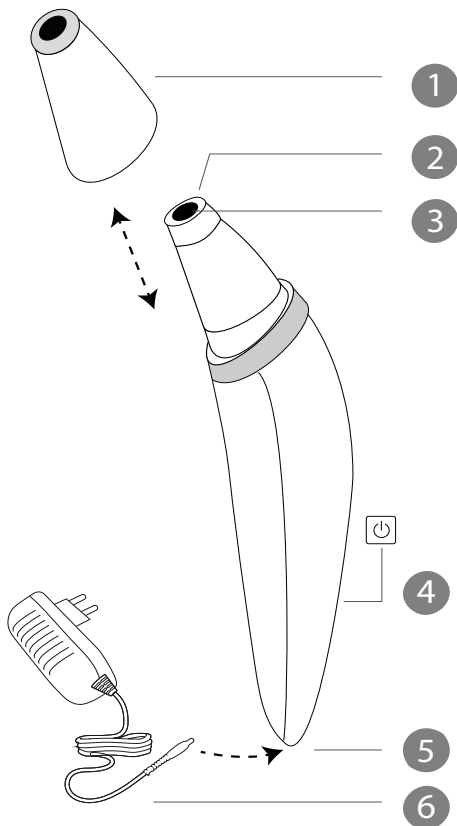Legal Notice

Home Skinovations Ltd. reserves the right to make changes to its products or specifications to improve performance, reliability, or manufacturability. Information furnished by Home Skinovations Ltd. is believed to be accurate and reliable at the time of publication. However, Home Skinovations Ltd. assumes no responsibility for its use. No license is granted by its implication or otherwise under any patent or patent rights of Home Skinovations Ltd.

No part of this document may be reproduced or transmitted in any form or by any means, electronic or mechanical, for any purpose, without the express written permission of Home Skinovations Ltd. Data is subject to change without notification.

Home Skinovations Ltd. has patents and pending patent applications, trademarks, copyrights, or other intellectual property rights covering subject matter in this document. The furnishing of this document does not give you any license to these patents, trademarks, copyrights, or other intellectual property rights except as expressly provided in any written agreement from Home Skinovations Ltd. Specifications are subject to change without notice.

SILK'N and the SILK'N logo are registered trademarks of Home Skinovations Ltd. Tavor Building, Shaar Yokneam POB 533, Yokneam 2069206, ISRAEL, www.silkn.eu or www.silkn.com, info@silkn.com.

Table of Contents

1.	Get to know your Silk'n ReVit Essential.....	19
2.	Contraindications & warnings	20
3.	Do not use Silk'n ReVit Essential if:	22
4.	What to expect with Silk'n ReVit Essential.....	23
5.	Possible side effects.....	23
6.	Before your first treatment	24
7.	Get started!.....	25
8.	Maintenance and storage.....	26
9.	Spare parts and consumables.....	27
10.	Troubleshooting	28
11.	Customer service.....	28
12.	Warranty	29
13.	Disposal	29
14.	Specifications	29
15.	Signs and symbols	30

1. Get to know your Silk'n ReVit Essential

Silk'n ReVit Essential is an electrical device. As such, it should be used with special attention to safety. Before using the device for the first time, be sure to read all warnings, contraindications and safety precautions. It is important that you read and understand the instructions on how to use the device including after-use procedures, and follow them strictly. We recommend you refamiliarise yourself with the instructions before treatment.

1.1 Package contents

Your Silk'n ReVit Essential package includes these items.

- ReVit Essential device
- Pack of filters
- 2 Treatment tips (Fine and Coarse)
- Spare rubber caps
- AC/DC adapter with cable
- Storage bag
- User manual
- Warranty card

1.2 Device description

ReVit Essential is a dual-action therapy device designed for home use. It combines microdermabrasion exfoliation and vacuum stimulation of the skin.

Microdermabrasion exfoliation

- Gentle, non-abrasive exfoliation to allow for smoother, fresher skin.
- Helps improve skin tone, texture and elasticity.
- Helps reduce the appearance of age spots.
- Helps reduce the appearance of fine lines around the eyes.

Vacuum stimulation

- Helps stimulate for smoother, more youthful looking skin.
- Brings some impurities to the surface.
- Collects the dead skin cells for disposal.

The device has a *treatment tip* (#1), *rubber cap* (#2), *filter*, (#3), *control button* (#4), *device socket* (#5), and an *adapter with cable*, (#6). Refer to the diagram at the front of this manual to see the location of each feature on the device.

1.3 Intended use

EN

Silk'n ReVit Essential is an over the counter device intended for the exfoliation of the skin's upper layer, for a fresher appearance. It is intended for the treatment of:

- Ageing skin of the face, neck and hands
- Fine facial lines
- Skin discolouration
- Dark spots and shallow scars
- Dull and sun-damaged skin

The device is intended for self-treatment at home. Do not use for any purpose other than as described in this manual.

2. Contraindications & warnings

Be sure to read this section carefully and take note of the warnings and contraindications.

- ! This device must not be used by: children, people under the age of 16, those with limited physical, sensory or psychological capacities, or people lacking experience or knowledge in how to use the device. Children should be supervised to ensure that they do not play with the device.

Do not use Silk'n ReVit Essential if you have:

- ! Rosacea.
- ! Poorly controlled endocrine disorders, such as diabetes and thyroid conditions.

- ❗ Abnormal skin conditions caused by diabetes, systemic or other metabolic diseases.
- ❗ Received radiation therapy or chemotherapy treatments within the past 3 months, or have areas of potential skin malignancies.
- ❗ Used physician-prescribed drugs on your skin within the past 6 months.
- ❗ Any other condition which, in the opinion of your physician or dermatologist, would make it unsafe for you to be treated.

Do not treat areas that have:

- ❗ Active eczema, psoriasis, lesions, open wounds or infections (cold sores), moles or warts, abnormal skin conditions caused by systemic or metabolic diseases (diabetes for example).
- ❗ A history of herpes outbreaks.
- ❗ Traumatized or injured skin.
- ☀️ Had recent exposure to the sun without protection (wait 1 day); for tanned or sunburnt skin, wait 3 days.

Additional warnings

- ❗ Do not use on nipples or genitals.
- ❗ Recommended use is no more than twice a week.
- ❗ Use the low suction mode if you experience any discomfort.
- ❗ Stop use immediately if your skin is very red.
- ❗ Remove all jewellery before treatments.

-  Protect your skin with SPF 30 or higher if you expose your skin to the sun after a treatment.
-  Always ensure that your hands and skin are dry before operating the device or touching the adapter plug.
-  Keep your device away from water or liquids! Do not use near bathtubs, showers, basins or other vessels containing water. Do not place or store the device where it can fall, be pushed or placed anywhere near water or any other liquid. There is a risk of electrocution. Unplug immediately and do not reuse.
-  Never leave the device unattended when plugged into an outlet. Keep the AC/DC adapter accessible at all times, so it can be removed from the socket easily in case of emergency. Unplug the device immediately after use.
-  Do not carry the device by the AC/DC adapter cable.
-  If the device malfunctions, unplug immediately.
-  Do not attempt to modify, open or repair your Silk'n ReVit Essential device. You may expose yourself to dangerous electrical components and may cause serious harm to yourself.

3. Do not use Silk'n ReVit Essential if:

- The adapter, attachments or accessories are not recommended by Home Skinovations Ltd.
- The device has a damaged cable or plug. Keep the power cable away from heated surfaces.
- You see or smell smoke when it is in use.
- The device is not working properly, or it appears to be damaged.
- The device casing is cracked or coming apart.

- The device has fallen into water.

In such cases, contact Silk'n Customer Service.

4. What to expect with Silk'n ReVit Essential

During a treatment session it is normal to experience or feel:

- **A humming sound from the device:** This is the suction pump working.
- **Some slight redness of the skin:** This should subside within a few hours, or up to one day.
- ❗ Stop use immediately if your skin is very red.
- ❗ Contact your physician, if you experience any complications.

5. Possible side effects

If you follow the instructions fully and pay attention to all warnings and contraindications, complications and side effects are uncommon. All cosmetic procedures, including those designed for home use, involve some degree of risk. Therefore, it is important that you understand and accept the risks and complications that may occur.

Adverse reaction	Risk
Skin redness You may experience dilation of the superficial blood capillaries. This may result in some redness for a few days. Additionally, some capillaries may present red spots that will also subside with time. In such cases, immediately cool the area with cold water.	Minor
Bruising You may experience slight bruising which should disappear after 2-3 days.	Minor

Consult your physician if you experience any persistent symptoms.

6. Before your first treatment


To use your Silk'n ReVit Essential, follow these steps:

1. Remove the device and other components from the box.
2. Plug the *adapter cable* into the device socket and plug the *adapter* into an electrical outlet.
3. Review the recommendations on how to use ReVit Essential on different treatment areas.

6.1 Treatment areas

We recommend that you start with the low suction mode and move to the high mode after a few uses. If you have any discomfort, work the treatment area in low suction mode.

Treatment area	Recommendations
Forehead, cheeks and chin	Move the device over your skin from the centre outwards and towards the hairline.
Eyes and nose	Use short upward strokes for the area around your eyebrows. Avoid the eyelids. For the nose area, use straight downward strokes following the curve of the nose and the contour of the cheekbones.
Neck	Use smooth, downward strokes and work the whole of the neck area.

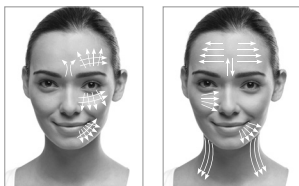To start with, you may work each area 2-4 times. When your skin is used to the treatment, you may increase to 4-6 passes.


6.2 ReVit Essential's Fine and Coarse tips

The device has two treatment *tips*. We recommend that you start with the Fine tip.

- Use the Fine *tip* for general treatment or delicate finishing.
- Use the Coarse *tip* for heavy duty exfoliation.

To change tips:

1. Make sure the device is not connected to an electrical outlet.
2. Remove the current *tip* (#1) from the device and replace with the other *tip*. Make sure the *tip* is firmly in place before use.

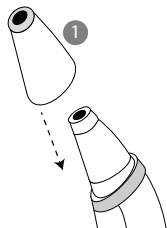7. Get started!

You're on your way to a fresher, smoother skin! Treating with Silk'n ReVit Essential is easy. Just follow these steps for each area you wish to treat. A total treatment should take around 5 minutes.

7.1 Recommended precautions

- Before you treat an area for the first time, do a spot test on a small area of skin. If you have no adverse reaction, go ahead and treat that area.
- Use the low suction mode for your first 2-3 treatments.
- We recommend you do not use the device more than twice a week.
- ! Stop immediately if you feel any discomfort.

7.2 Treatment steps

1. Make sure your skin is clean, dry and free from powder, cream or gel.
2. Make sure that the *rubber cap* (#2), *filter* (#3), and *tip* (#1) are firmly in place. For your first few treatments, be sure to use the *Fine tip*.
3. Press the *control button* (#4) to turn on the device. Your ReVit Essential is ready for use in low suction mode.
4. Apply the *tip* firmly to your skin so that you feel the suction.
5. Work each area as recommended in section "6.1 Treatment areas".
6. To move from low to high suction mode, press the *control button* again.
7. To switch off the device after treatment, short press the *control button* once or twice, depending on which suction mode you are in.
8. Unplug the device and store ReVit Essential in its original box. Please refer to section "8. Maintenance and storage".

7.3 After treatment

- Avoid direct sunlight immediately after treatment, and apply sun factor SPF 30 or higher.
- You may wish to apply a facial mask, or apply a day or night cream.

8. Maintenance and storage

It is important to clean your Silk'n ReVit Essential device and store it properly after each treatment.

First, unplug the device.

To clean the device, wipe with a clean, dry cloth. Use a cleaning fluid with no more than 70% alcohol, or another non abrasive cleaning fluid. Never use acetone or other abrasive fluids, because they will damage your device.

To clean the *tip*, first remove it from the device, wash with soap and water, then dry thoroughly.

To clean the *filter*, first remove it from the device, brush away any debris using a soft brush, and place it back on the device.

After cleaning, we recommend you store the device in its original box.


With the exception of the *tips*, never immerse your Silk'n ReVIt Essential in water.

8.1 Replacing the filter

The *filter* needs cleaning or replacing due to the build-up of discarded skin. We recommend that you do this every few treatments.

To replace the filter:

1. Remove the *tip* from the device . You will then see the *rubber cap* and the small black *filter*.
2. Remove the *rubber cap* (#2), and then remove the *filter* (#3).
3. To clean the *filter*, brush off any debris, then place it back on the device. To change a *filter*, remove the current *filter* and insert a new one.
4. Place the *rubber cap* back into position, then put the *tip* back on the device, making sure that it is firmly connected.

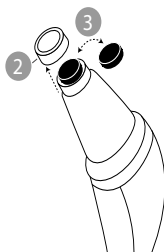8.2 Replacing the tip

We recommend that you replace the *tip* every six months, or earlier if the abrasive end of the tip is damaged.

9. Spare parts and consumables

You can order more *filters* and *tips* from our website: www.silkn.eu. For replacement *rubber caps* or *AC/DC adapters*, please contact Silk'n Customer Service.

10. Troubleshooting

Do not attempt to open or repair your Silk'n ReVit Essential. You may expose yourself to dangerous electrical components and cause serious harm to yourself. It will also void the warranty.

Problem	Checks
My Silk'n ReVit Essential does not start.	Make sure the <i>adapter</i> is properly connected to the device, that it is plugged into an electrical outlet and that you have pressed the <i>control button</i> .
I don't feel any suction when the device touches my skin.	Make sure the device is on, and check that the <i>filter</i> , <i>rubber cap</i> and <i>tip</i> are placed on the device correctly.

Only authorized personnel are permitted to perform repairs. For problems or queries, please contact Silk'n Customer Service.

If the device fails, unplug immediately.

11. Customer service

See our website for more information: www.silkn.eu. This manual is also available as a PDF download from: www.silkn.eu.

Contact Silk'n Customer Service to report unexpected operations or events, or any other problem with your device. Our contact details are:

Country	Service Number	Email
Belgium	0900-25006	serviceconsommateurbe@silkn.eu
France	0891-655557	serviceconsommateurfr@silkn.eu
Germany	089 51 23 44 23	kundenservice@inno-essentials.de
Italy		servizioclientiit@silkn.eu
Netherlands	0900-2502217	klantenservicenl@silkn.eu
Spain	900 823 302	servicioalcliente@silkn.eu

Country	Service Number	Email
UK	0906-2130009	customercareuk@silkn.eu
US and Canada	1-877-367-4556 (toll free)	contact@silkn.com
Other countries	+31 (0)180-330 550	info@silkn.com

12. Warranty

Refer to your Silk'n warranty card.

13. Disposal

Help protect the environment. Do not dispose with household waste and be sure to follow your country's laws for the disposal of electrical items.

14. Specifications

Model no.	H4201
Technology	Air vacuum filling tip
Input power	12 V DC, 0.2A
System weight	110 gm
Transport & storage between uses and storage condition	Temperature: -40°C to +70°C Relative humidity : 10% to 90% Atmospheric pressure: 500 to 1060 hPa
Operating conditions	Temperature: 10°C to 35°C Relative humidity: 30% to 75% Atmospheric pressure: 700 to 1060 hPa

Adapter models and rating	SAW12-120-0500XX XX Refer to plug type Input: 100-240VAC, 50-60 Hz, 0.3A Output: 12V, 0.5A
Expected service life of the equipment	3 years
Compliance	2014/35/EU 2014/30/EU Directive 2011/65/EC RoHS REACH regulation EC No. 1907/2006

15. Signs and symbols

	Products marked with this symbol fulfil the requirements of the EU Directives.
	Class III. The appliance is supplied from an extra-low voltage power source.
	For indoor use only.
	Do not use in wet environments.
	Do not dispose with household waste. Dispose in accordance with your country's laws and help protect the environment. In accordance with the WEEE Directive.

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InnoEssentials Deutschland GmbH